

SOLOCARE AQUA™




EASY CLEAN

Comfort and Convenience

Free MicroBlock™ Antibacterial Lens Case Included

- Clinically proven to reduce lens case contamination¹
- Helps to kill bacteria
- Resists the growth of new bacteria on the lens case

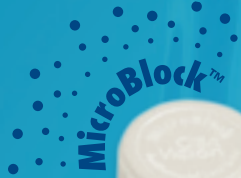

Together, Sorbitol and Dexpantenol Lock in Moisture and Help Prevent Lens Dehydration

- In a comparative in vitro study, lenses soaked in SOLOCARE AQUA retained up to 85% more moisture than lenses soaked in ReNu® MultiPlus®, COMPLETE® and OPTI-FREE® EXPRESS®²

With SOLOCARE AQUA, you can help patients maximize the silicone hydrogel contact lens-wearing experience.

SOLOCARE AQUA is the first multi-purpose solution with the exclusive MicroBlock lens case, to reduce bacterial contamination of the lens case. The plastic of the lens case is infused with silver ions—an antibacterial agent—to destroy many strains of bacteria and help prevent growth of bacteria in the lens case.

SOLOCARE AQUA also offers a convenient 5-minute option to clean lenses easily whenever needed.*

Ideal for soft contact lens wearers who:

- Lead an active life
- Need a quick clean and don't have the 4 to 6 hours required by other solutions
- Sleep in their lenses**
- Wear their lenses for extended periods of time during the day

* With rub. ** Only an eye care professional can determine if a patient may wear their prescribed lenses overnight.

From the Worldwide Leader
in Silicone Hydrogel Technology

**CIBA
Vision**
A Novartis Company

What's Inside Makes the Difference

INGREDIENTS	SOLOCARE AQUA
Disinfectant (Active Ingredient) Kills harmful microorganisms	Polyhexanide (0.0001%)
Buffer/Comfort Maintains pH Balances solution for gentleness and comfort	Tromethamine Sodium Phosphate
Surfactant/Cleaning Loosens dirt, debris and deposits	Poloxamer 407
Moisturizing Agent(s) Attracts Moisture	HydroLock Effect Dexpanthenol Sorbitol

HydroLock Formulation

As a moisturizer, dexpanthenol creates a barrier effect which helps retard the evaporation of moisture. Sorbitol is a humectant—a substance that absorbs and promotes the retention of moisture.



Directions For Use

SOLOCARE AQUA can be used in **2 Convenient Ways**

1 No Rub Directions for Use:

- Fill the lens case with enough solution to cover the lenses. Store the lenses in the closed lens case for 4 hours to overnight

2 5-Minute Express Option (with a Rub):

- Apply at least 3 drops of SOLOCARE AQUA on each side of the lens and gently rub the lens for 10 seconds per lens side
- Rinse the lens thoroughly with SOLOCARE AQUA
- Fill the lens case with fresh solution. Allow your lenses to soak for at least 5 minutes



Offer Your Silicone Hydrogel Patients

COMFORT AND CONVENIENCE



ReNu, MoistureLoc and MultiPlus are registered trademarks of Bausch & Lomb, Inc. OPTI-FREE and EXPRESS are registered trademarks of Alcon Laboratories, Inc. COMPLETE is a registered trademark of Advanced Medical Optics, Inc.

References: 1. CIBA Vision. Data on file, 2005. 2. CIBA Vision. Data on file, 2004.

© 2006 CIBA Vision Corporation 2006-03-0180

CIBA Vision
A Novartis Company